

ASSEMBLY INSTRUCTION
OAO-86A HCP Fixing Bracket for DS-8100



FUKUDA DENSHI CO., LTD.

3-39-4, Hongo, Bunkyo-ku, Tokyo, Japan
Tel: +81-3-5684-1455 Fax: +81-3-3814-1222

Printed in Japan 4L0110690 201403

Copyright © 2014 by Fukuda Denshi Co., Ltd.

No part of this document may be copied or transmitted in any form without the prior written permission of Fukuda Denshi Co., Ltd.

- Thank you for purchasing this product.
- Before using this product, read through this assembly instruction and also refer to the DS-8100 operation manual.
- After reading, keep this assembly instruction near the unit for future reference.



This device bears the CE label in accordance with the provisions of Medical Device Directive 93/42/EEC and RoHS Directive 2011/65/EU.



Fukuda Denshi UK
Unit 7, Genesis Business Park, Albert Drive, Woking, Surrey GU21 5RW, United Kingdom

This product is used for fixing the HCP-810 CO₂ Gas Unit or the HPD-810 Gas Unit I/F to the DS-8100 monitor. Follow this assembly instruction, and use this product correctly.

Safety Precautions

Make sure to follow the precautions indicated below, as these are important messages related to safety. The warning displayed and symbols used in this instruction are shown below.

	CAUTION	Failure to follow this message may cause injury or failure of the product.
NOTE		A note is not related to product safety, but provides information about the correct use and operating procedures to prevent incorrect operation and malfunction of the equipment.

NOTE


For quality improvement, specifications are subject to change without prior notice.

Contents

Refer to the following list and verify that all items are included.

HCP Holder x1	HCP Bracket x1	Screw M3X8 x4

Example

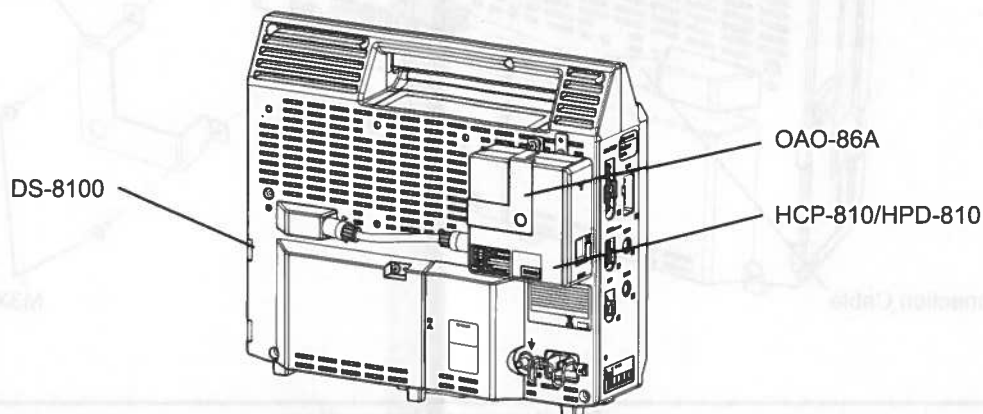

CAUTION

- This product is used when fixing the HCP-810/HPD-810 to the DS-8100 series patient monitor. Do not use this product for any other purpose.
- Check that all the screws are tightened correctly after assembly and at periodic inspections.

Assembly Procedure


This procedure explains how to attach OAO-86A to the monitor. The process requires the following tool.

- Phillips screwdriver #2


NOTE

- Attach this product to the monitor on a level surface following the procedure below.
- A Phillips screwdriver is required to complete the installation.
- Make sure that the screws are securely tightened.

1. Attach the HCP holder to the back of the DS-8100 main unit using the screws.

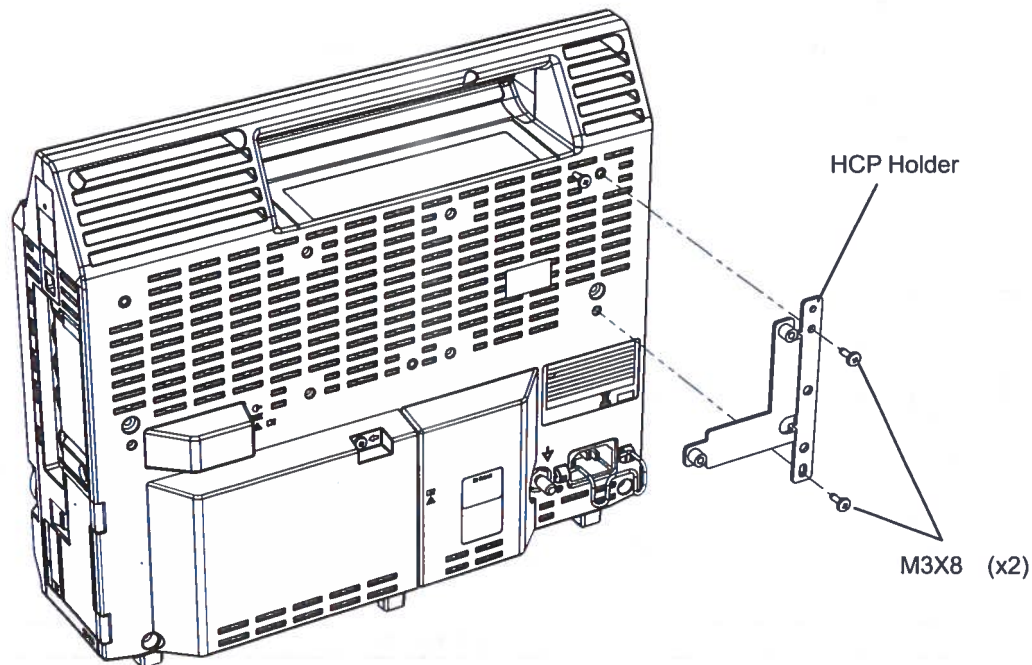

2. Place the HCP-810/HPD-810 on the HCP holder and attach the HCP bracket using the screws.

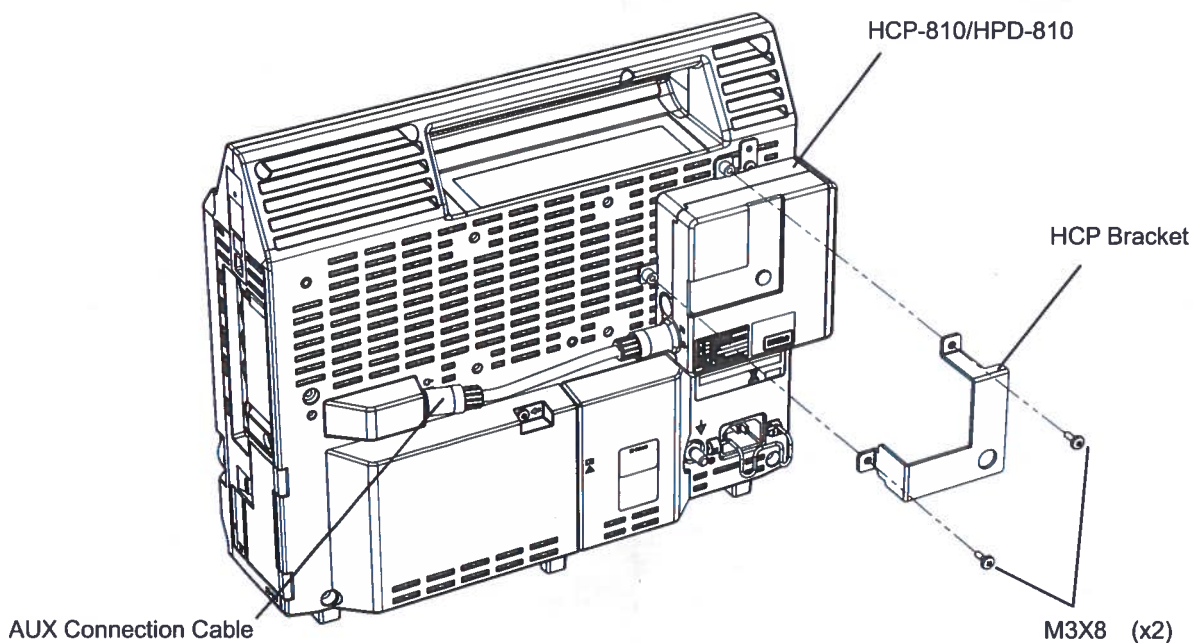

Maintenance / Inspection

Verify that the screws are securely tightened once a year.