

# Operation Manual OA-470 Relay Cable Mounting Bracket

**FUKUDA DENSHI CO., LTD.**

39-4, Hongo, 3-chome, Bunkyo-ku, Tokyo, Japan  
Phone: +81-3-3815-2121 Fax: +81-3-3814-1222

Printed in Japan 4L3537 200502



Copyright © 2005 by Fukuda Denshi Co., Ltd.  
No part of this document may be copied or transmitted in any form without the prior written permission of Fukuda Denshi Co., Ltd.

- Thank you for purchasing this product.
- Before using, read this operation manual thoroughly.
- Keep this operation manual near the product where it can be always referred.

**CE** This device bears the CE label in accordance with the provisions of  
**Medical Device Directive 93/42/EEC.**

THE PERSONS RESPONSIBLE FOR PLACING DEVICES ON THE EC MARKET UNDER MDD 93/42/EEC  
NAME : FUKUDA DENSHI UK  
ADDRESS : 13 WESTMINSTER COURT, HIPLEY STREET OLD WOKING, SURREY GU22 9LG, U.K.

This product is a mounting bracket to attach the relay cable to the HS-700 Super Module.  
Follow the procedure explained in this operation manual for proper attachment.

## Safety Precautions

Read the following precautions thoroughly to ensure proper, correct and safe use of the product.  
The following is a description of the safety and precaution messages used in this manual.

<b>CAUTION</b>	Ignoring these instructions may cause injury or damage to the equipment.
<b>NOTE</b>	A note is not related to product safety. A note provides information about the proper use and operating methods to prevent incorrect operation and trouble with the unit.

## Package Composition

Refer to the following list and verify the item. If any shortage or damage is found, please contact Fukuda Denshi.


	Item	Q'ty
1	OA-470 Relay Cable Mounting Bracket	1
2	Screw M3 × 8	2

## Procedure


### NOTE

- Attach this product on a level surface according to the procedure explained below.
- Prepare a Phillips screwdriver for fixing on the OA-470.

1. Disconnect the 2 screws attached to the panel of the HS-700. Fix the OA-470 (1) to the HS-700 using the accessory screws (2). (It can be fixed on to the opposite side using the same procedure.)  
The detached screws can be discarded as these will not be necessary.

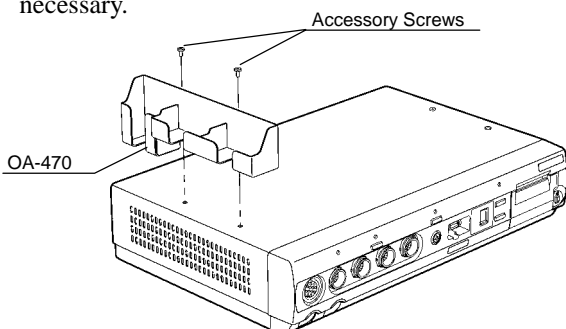


2. If using in conjunction with the OA-469 (DS-7300 Mounting Bracket), fix on with the same procedure as above using the accessory screws (2).

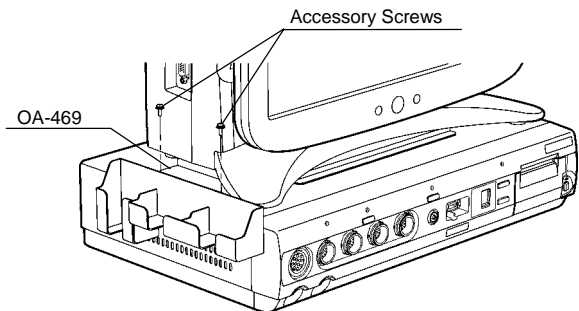

### 【Example of Usage】

The 2ch BP relay cable (CJ-7546) and 2ch temperature relay cable (CJ-7414) can be placed on the OA-470 as shown on the following illustration.

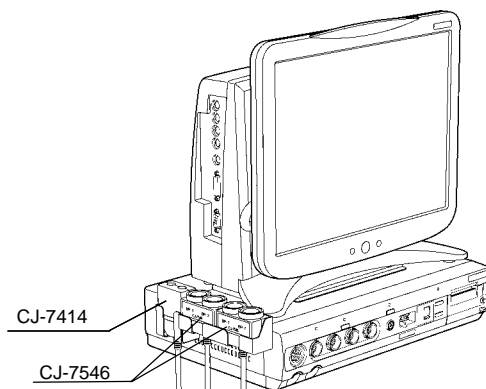

### CAUTION

- This product is an exclusive bracket for attaching the relay cable to the HS-700 Super Module. Do not use for other purposes.
- Verify secure connection of the screws after the attachment and at periodic inspection.