

Operation Manual OA-473 GCX Plate Adapter

FUKUDA DENSHI CO.,LTD.

39-4, Hongo, 3-chome, Bunkyo-ku, Tokyo, Japan
Phone:+81-3-3815-2121 Fax:+81-3-3814-1222



Printed in Japan 4L3599A 200502

Copyright © 2005 by Fukuda Denshi Co., Ltd.

No part of this document may be copied or transmitted in any form without the prior written permission of Fukuda Denshi Co., Ltd.

- Thank you for purchasing this product.
- Before using, read this operation manual thoroughly.
- Keep this operation manual near the product where it can be always referred.



This device bears the CE label in accordance with the provisions of Medical Device Directive 93/42/EEC.

THE PERSONS RESPONSIBLE FOR PLACING DEVICES ON THE EC MARKET UNDER MDD 93/42/EEC

NAME : FUKUDA DENSHI UK

ADDRESS : 13 WESTMINSTER COURT, HIPLEY STREET OLD WOKING, SURREY GU22 9LG, U.K.

This product is an adapter used for attaching the GCX plate to the HS-700 series.

Safety Precautions

Read the following precautions thoroughly to ensure proper, correct and safe use of the product.

The following is a description of the safety and precaution messages used in this manual.

CAUTION	Ignoring these instructions may cause injury or damage to the equipment.
NOTE	A note is not related to product safety. A note provides information about the proper use and operating methods to prevent incorrect operation and trouble with the unit.

Package Composition

Refer to the following list and verify the item. If any shortage or damage is found, please contact Fukuda Denshi.


No.	Item	Q'ty
1	Adapter	1
2	Screws (M3x8)	4

Procedure

NOTE


- The attachment of this product should be performed on a level surface according to the procedure explained below.
- Prepare a Phillips screwdriver for fixing on the adapter.

1. Detach the 4 rubber supports from the bottom of the equipment. The detached rubber supports and screws can be discarded as these will not be necessary.

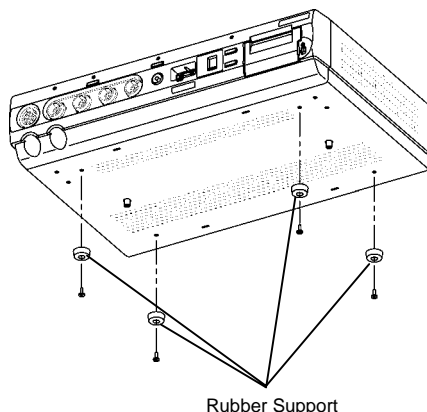


Rubber Support

2. Fix on the adapter (1) to the bottom of HS-700 using the accessory screws (2).

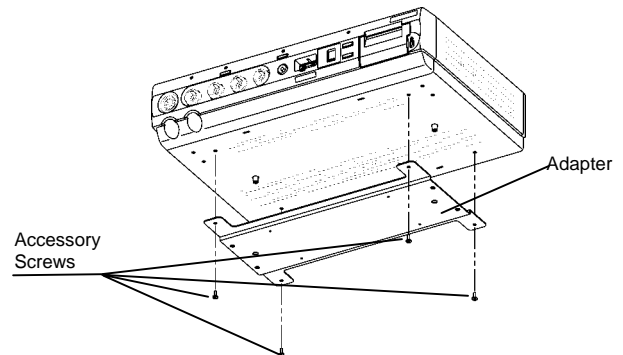


【Final Drawing】

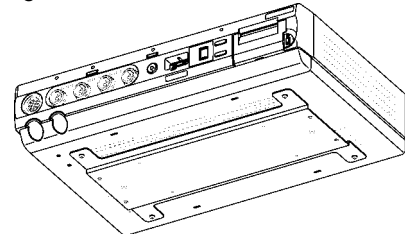

【Reference】GCX Plate Attachment

Fix on the GCX plate to the bottom of the equipment using 4 screws (accessory of GCX plate).

The GCX plate is not included in the package. Purchase it separately.

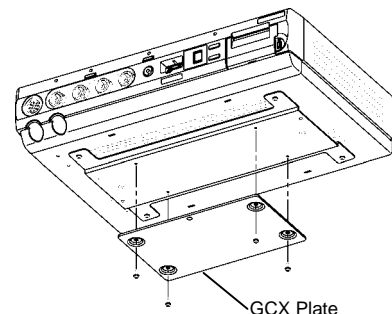

GCX Plate

CAUTION

- This product is exclusively intended for attaching the GCX plate to the HS-700 series. Do not use for other purpose.
- If attaching the equipment to the wall mount, check the specification of the wall mount (withstand load, etc.) in advance. Also, consult with the construction company about the wall material, thickness, etc.
- Verify secure connection of the screws after the attachment and at periodic inspection.