

# Operation Manual OA-469 DS-7300 Mounting Bracket

**FUKUDA DENSHI CO., LTD.**

39-4, Hongo, 3-chome, Bunkyo-ku, Tokyo, Japan  
Phone: +81-3-3815-2121 Fax: +81-3-3814-1222

Printed in Japan 4L3536A 200502



Copyright © 2005 by Fukuda Denshi Co., Ltd.  
No part of this document may be copied or transmitted in any form without the prior written permission of Fukuda Denshi Co., Ltd.

- Thank you for purchasing this product.
- Before using, read this operation manual thoroughly.
- Keep this operation manual near the product where it can be always referred.



**This device bears the CE label in accordance with the provisions of Medical Device Directive 93/42/EEC.**

THE PERSONS RESPONSIBLE FOR PLACING DEVICES ON THE EC MARKET UNDER MDD 93/42/EEC

NAME : FUKUDA DENSHI UK

ADDRESS : 13 WESTMINSTER COURT, HIPLEY STREET OLD WOKING, SURREY GU22 9LG, U.K.

This product is a mounting bracket to attach the DS-7300 Bedside Monitor to the HS-700 Super Module.  
Follow the procedure explained in this operation manual for proper attachment.

## Safety Precautions

Read the following precautions thoroughly to ensure proper, correct and safe use of the product.

The following is a description of the safety and precaution messages used in this manual.

### CAUTION

Ignoring these instructions may cause injury or damage to the equipment.

### NOTE

A note is not related to product safety.  
A note provides information about the proper use and operating methods to prevent incorrect operation and trouble with the unit.

## Package Composition

Refer to the following list and verify the item. If any shortage or damage is found, please contact Fukuda Denshi.


	Item	Q'ty
1	OA-469 Mounting Bracket	1
2	Screw (M3 × 8)	8
3	Inch Screw	4

## Procedure


### NOTE

- Attach this product on a level surface according to the procedure explained below.
- Prepare a Phillips screwdriver for fixing on the mounting bracket.

1. Fix on the OA-469 Mounting Bracket (1) to the bottom of the DS-7300 using the accessory screws (2).

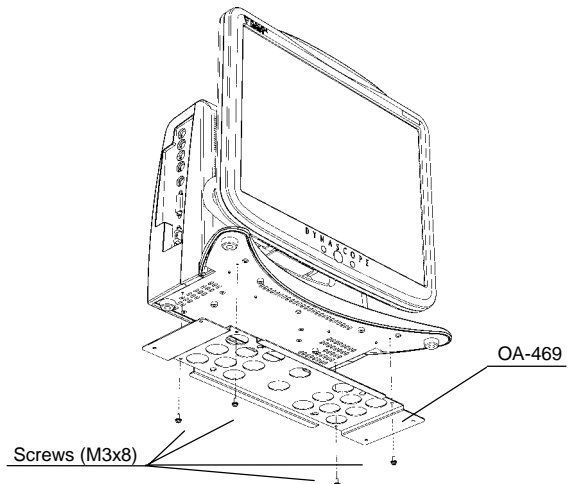

2. Detach the 4 screws from the top panel of the HS-700. Place the OA-469 attached DS-7300 on the HS-700, and fix on with the accessory screws (2). The detached screws can be discarded as these will not be necessary.

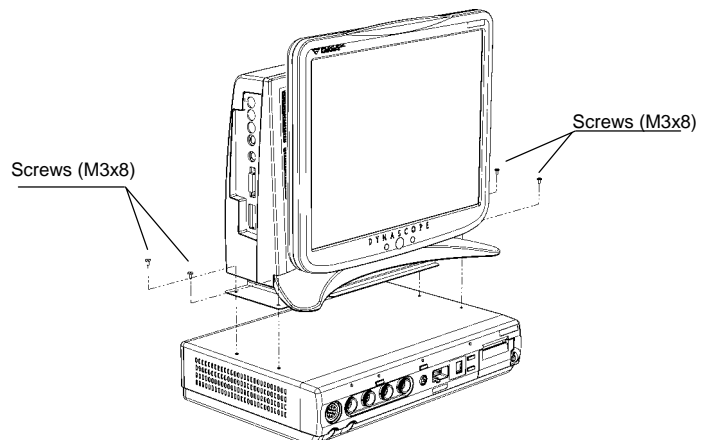

### CAUTION

- This product is exclusively intended for attaching the DS-7300 Bedside Monitor to the HS-700 Super Module. Do not use for other purposes.
- Verify secure connection of the screws after the attachment and at periodic inspection.