

# ASSEMBLY INSTRUCTION

## OA0-66A Stand for DS-8100



### FUKUDA DENSHI CO.,LTD.

3-39-4, Hongo, Bunkyo-ku, Tokyo, Japan  
Tel:+81-3-5684-1455 Fax:+81-3-3814-1222

Printed in Japan 4L011066A 201403

Copyright © 2014 by Fukuda Denshi Co., Ltd.

No part of this document may be copied or transmitted in any form without the prior written permission of Fukuda Denshi Co., Ltd

- Thank you for purchasing this product.
- Before using this product, read through this assembly instruction and also refer to the DS-8100 operation manual.
- After reading, keep this assembly instruction near the unit for future reference.



This device bears the CE label in accordance with the provisions of Medical Device Directive 93/42/EEC and RoHS Directive 2011/65/EU.



Fukuda Denshi UK  
Unit 7, Genesis Business Park, Albert Drive, Woking, Surrey GU21 5RW, United Kingdom

This product is used for fixing the DS-8100 monitor to:

- (1) Shelf (Option OA0-87A is required)
- (2) GCX Stand

Follow this assembly instruction, and use this product correctly.

## Safety Precautions

Make sure to follow the precautions indicated below, as these are important messages related to safety. The warning displayed and symbols used in this instruction are shown below.


	<b>CAUTION</b>	Failure to follow this message may cause injury or failure of the product.
<b>NOTE</b>		A note is not related to product safety, but provides information about the correct use and operating procedures to prevent incorrect operation and malfunction of the equipment.

## NOTE


For quality improvement, specifications are subject to change without prior notice.

## Contents

Refer to the following list and verify that all items are included.

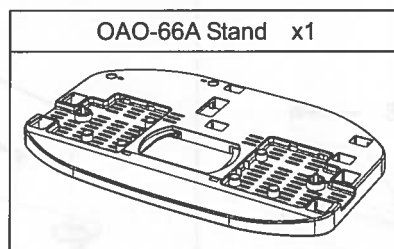

## Example

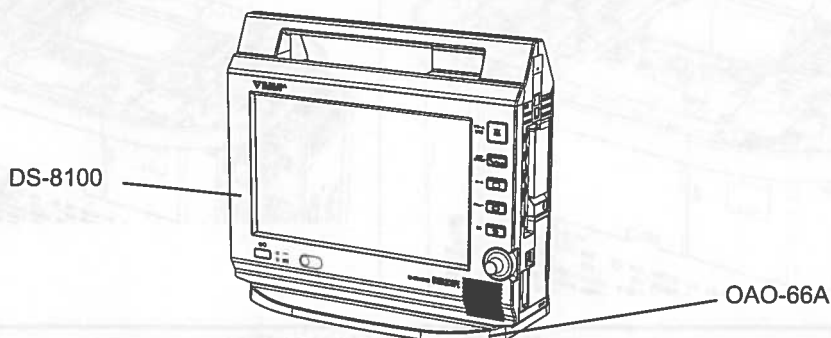

## CAUTION

- This product is used when fixing the DS-8100 series patient monitor to a shelf or the GCX stand. Do not use this product for any other purpose.
- Check that all the screws are tightened correctly after assembly and at periodic inspections.

## Assembly Procedure


This procedure explains how to attach OAO-66A to the monitor. The process requires the following tool.

- Phillips screwdriver #2


### NOTE

- Attach this product to the monitor on a level surface following the procedure below.
- A Phillips screwdriver is required to complete the installation.
- Make sure that the screws are securely tightened.


1. Remove the screw from the bottom side of the stand, then remove the stand bracket.

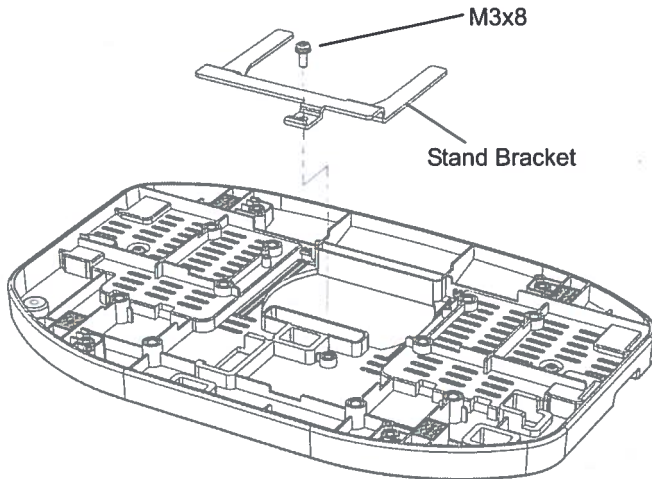

2. If the option units (HR-811, HR-810, CU-810) are used, attach them first on the back of the main unit.

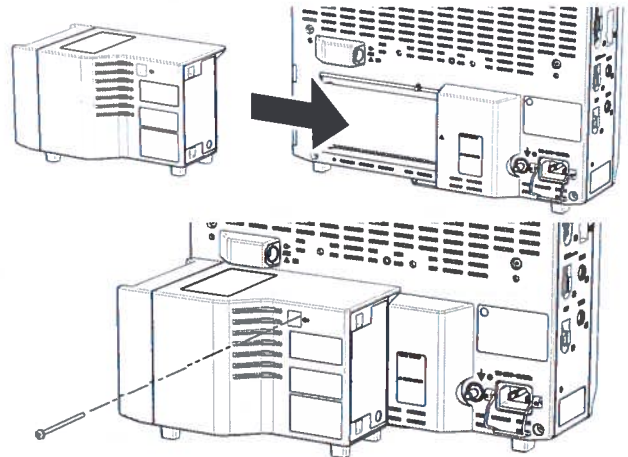


3. Place the stand on the bottom of the DS-8100 and attach the stand bracket using the screw removed in 1. Tighten the screws on each side.

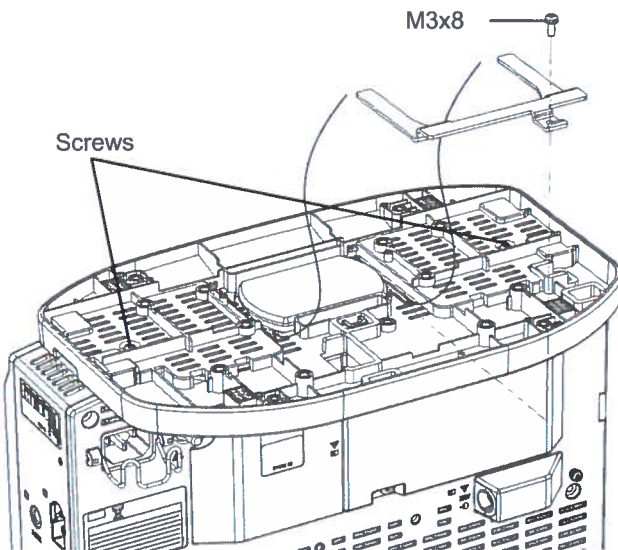

### <How to attach to GCX plate>

Secure the GCX plate (sold separately) with the 4 screws.

M4x12 (Double Washer Sems)

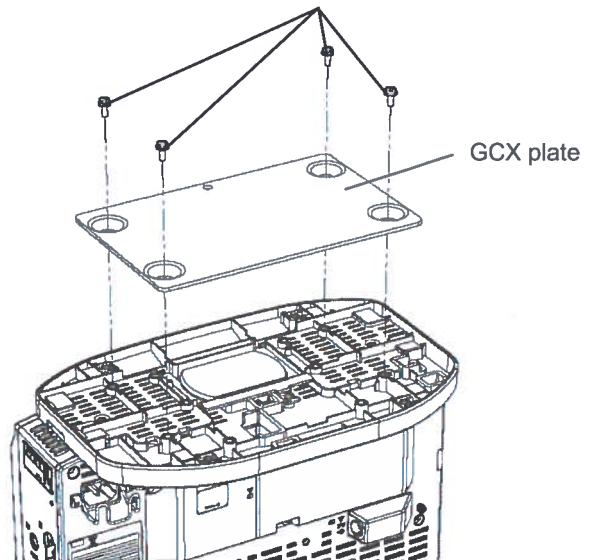

## Maintenance / Inspection

Verify that the screws are securely tightened once a year.