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Printed in Japan 4L010788B 201406

© 2014 Fukuda Denshi Co., Ltd.

No part of this document may be copied or transmitted in any form without the prior written permission of Fukuda Denshi Co., Ltd.

- Thank you for purchasing this product.
- Before using, read this instruction manual and the operation manual for the patient monitor thoroughly.
- Keep this manual near the product where it can be always referred to.

This device bears the CE label in accordance with the provisions of Medical Device Directive 93/42/EEC and RoHS Directive 2011/65/EU.

Fukuda Denshi UK  
Unit 7, Genesis Business Park, Albert Drive, Woking, Surrey GU21 5RW, United Kingdom

This product is a special output box which phone plugs can be inserted to extract the analog output from the HS-8000 Super Unit. By using this product, ECG (1ch), BP (2ch), and the QRS synchronized signal can be output.

## Safety Precautions

Make sure to follow the precautions indicated below, as these are important messages related to safety.

<b>⚠ WARNING</b>	Failure to follow this message may result in death or serious injury.
<b>⚠ CAUTION</b>	Failure to follow this message may cause injury or failure to the equipment.
<b>NOTE</b>	A note is not related to product safety, but provides information about the correct use and operating procedures to prevent incorrect operation and malfunction of the equipment.

## Contents

Refer to the following list and verify that all items are included.

	Item	Qty	Note
1	CJO-C01Q-SJ0.3	1	—
2	Clip	1	—
3	Instruction Manual	1	—


## Graphic Symbol

Refer to the following for the meaning of the symbol indicated on the equipment.


Symbol	Meaning
	<b>Caution:</b> refer to accompanying documents Indicates the need to refer to related accompanying documents before operation.
	<b>Signal Output</b> Indicates the connector which outputs signals (analog, and QRS sync.).

## Connection Procedure

<Connection to the HS-8000 Super Unit>
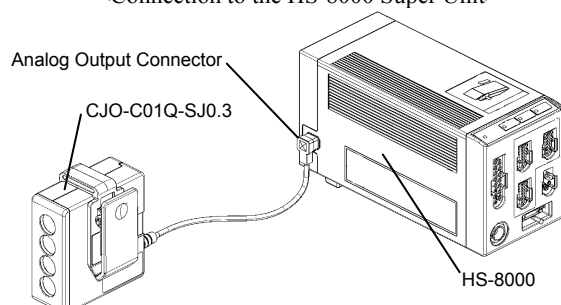
## Name of Parts and Their Functions

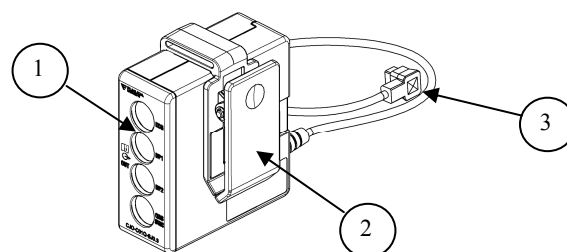	Name	Function
1	Analog Output Connectors	For inserting the phone plug to output the analog signal.
2	Clip	For attaching the output box to a bedside rail or headboard.
3	HS-8000 Connector	For connecting the output box to the HS-8000 series analog output connector.

### 1. Output Signals

ECG : ECG analog waveform (1 channel) can be output.  
BP1, BP2 : BP analog waveforms (2 channels) can be output.  
QRS SYNC : QRS synchronized signal can be output.

### 2. Analog Output Setup

For details of the analog output waveform, refer to the operation manual of the DS-8500 Patient Monitor and the HS-8000 Super Unit.

### ⚠ WARNING

- Only a phone plug should be connected to the connector. If any unintended cable (ex. temperature probe 400 series) is connected, a leakage current may cause electric shock to the patient and operator.
- Never remodel the equipment. It may cause a serious accident.

### ⚠ CAUTION

Do not input a signal to the connectors as there are all for output only. If done so by mistake, it may cause failure to the equipment.

### NOTE

For quality improvement, specifications are subject to change without prior notice.